
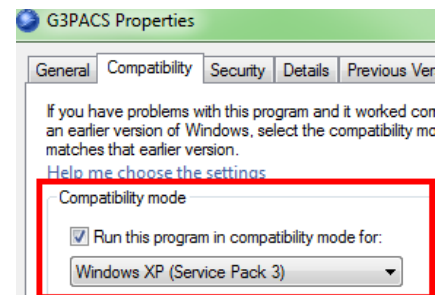


GOImaging - PACS Install Guide

<https://pih.infinittna.com>

You will need Administrator privileges on your computer to install INFINITT PACS.

- Open **Microsoft Internet Explorer** and enter the web address shown above.
 - Click the *Tools* menu and then choose *Internet Options*.
- On the **General** tab click on *Browsing History Settings*.
 - Choose *Every time I visit the webpage* and click *OK*.
- Go to the **Security** tab.
 - Select the *Trusted Sites* zone and then click on the *Sites* button.
 - Add this website to the zone: <https://pih.infinittna.com> - then click *Add*.
 - Click the *Custom level* button for the *Trusted Sites* zone.
 - Reset custom settings to **Low**. Hit the **Reset** button and then *OK* to exit.
- Go to the **Privacy** tab.
 - Turn OFF the pop-up blocker.
- **Internet Explorer 10 and higher:**
 - Click the Compatibility Mode button (broken page icon). 
- **Windows 7 and Vista Pro users:**
 - Go to the *Control Panel > User Accounts*. Change *User Account Controls* settings to Never Notify.
- **Windows 8 users:**
 - Go to the folder *C:\Infinitt\Viewer* and find the *G3PACS.exe* program.
 - Right click on the icon and select *Properties*.
 - Choose the *Compatibility* tab and select the box to *Run this program in compatibility mode for:*



You are ready to use INFINITT RIS/PACS. Enter your *User ID* and *Password* on the login screen. Smart Update will run and install the necessary software. A shortcut will be installed on the desktop for quick login access.

System Requirements:

Operating System:	<i>Windows XP Pro , Windows Vista Pro, Windows 7 (any version)</i>
Monitor resolution:	<i>1024 x 768 (minimum) , 1280 x 1024 and higher (recommended)</i>
Web Browser:	<i>Internet Explorer versions - 7 , 8, 9, 10 (Internet Explorer ONLY)</i>
Computer Specs:	<i>1GB RAM (minimum) , 10GB hard drive space , High Speed Internet</i>

